

**medela** 

mother's milk,  
everyday amazing™

# Medela breast pumps with Flex™ technology


A whole new pumping  
experience personalised to you



**NEW**

# Medela Flex™ technology breast pumps

A whole new pumping experience, personalised to you

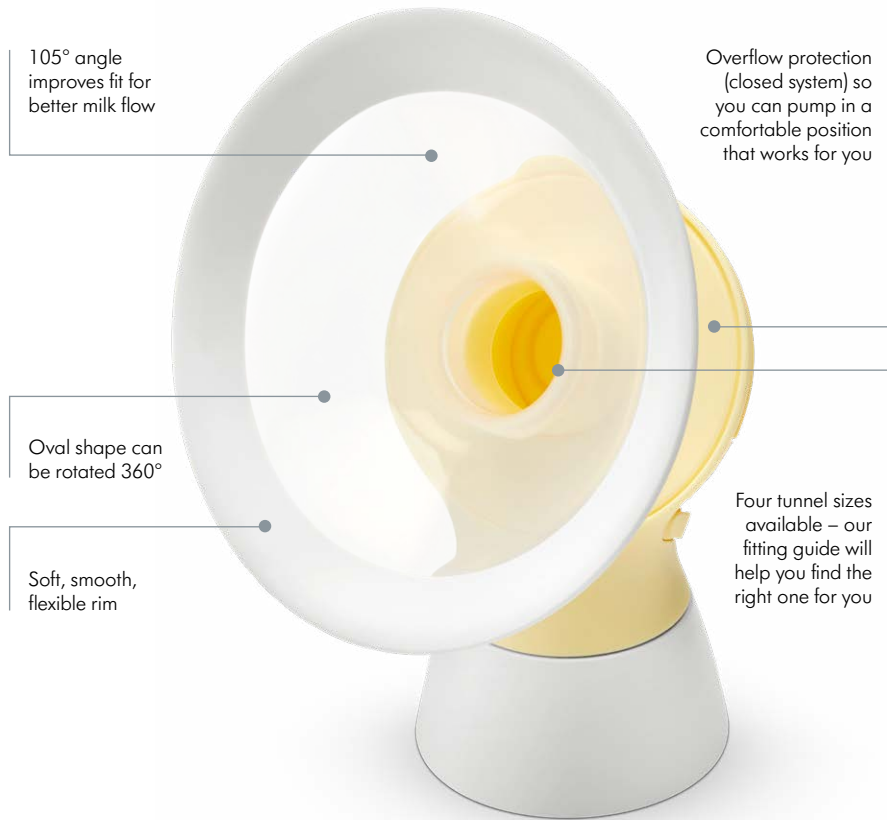

105° angle improves fit for better milk flow

Overflow protection (closed system) so you can pump in a comfortable position that works for you

Oval shape can be rotated 360°

Soft, smooth, flexible rim

Four tunnel sizes available – our fitting guide will help you find the right one for you

**Every mum is different, from the shape of her breasts and nipples to the demands of her lifestyle. Our new Flex™ technology breast pumps are the first designed to adapt to fit your shape for more efficient, comfortable expression.**



## The PersonalFit Flex™ breast shield

The new PersonalFit Flex™ breast shield was developed with the help of a team of research scientists and breastfeeding mums. It's designed to easily adapt to your shape, and encourages your milk flow by helping prevent compression of the milk ducts. This makes it easier for every mum to find the pumping position that's most comfortable for her. The result? More comfortable, more efficient expression, shaped around you.



## 2-Phase Expression® technology

When your baby feeds at the breast, he sucks quickly at first to start the milk flow, then more slowly to feed. Medela 2-Phase Expression® technology mimics these instinctive rhythms for a more natural-feeling and comfortable pumping experience.

## Double pumping

Research shows that the unique combination of double pumping with Medela 2-Phase Expression® technology delivers more milk, with a higher energy content, in less time than it takes to single express two bottles.

# More milk and comfort, as tested with mums

Our research scientists tested the PersonalFit Flex™ breast shield in four clinical studies with breastfeeding mums. Across more than 1,000 pumping sessions, they proved it delivered more comfort and efficiency. Like every Medela product, our Flex™ breast pumps are designed to make breastfeeding easier, so you can feed your baby with your amazing milk for longer.



**Mother's milk is amazing.** It's the best start in life your baby can have. Medela has been working with researchers in breast milk for more than 50 years. We're here to support every breast milk-feeding mother with our knowledge, expertise and unique range of products, including our Swiss-made breast pumps, the trusted choice of healthcare professionals worldwide.

# The Flex™ technology breast pump range

## Swing Flex™

For quick, efficient expression when you need it



## Swing Maxi Flex™

For that extra bottle in half the time



## Freestyle Flex™

The flexibility to pump any time, anywhere



Type of pump	Single electric	Double electric	Double electric
Recommended usage	Occasional	Daily	Daily
Easy-to-use and comfortable	•	•	•
PersonalFit Flex™ shield	•	•	•
Adjustable vacuum strength	•	•	•
Vacuum range	-45mmHg to -245mmHg	-45mmHg to -245mmHg	-45mmHg to -245mmHg
2-Phase Expression®	•	•	•
Over flow protection (closed system)	•	•	•
BPA free	•	•	•
Connects to MyMedela app			•
Backlit display showing session timer, battery status, vacuum level, bluetooth connection etc.			•
USB charger			•
Battery operated	Up to 1 hour (battery not included)	Up to 1 hour (battery not included)	Up to 8 double pumping session build-in rechargeable battery
All-in-one solution for storing, cooling and transporting expressed milk			Cooling bag, storage bottles and carry bag
Weight	200g w/o batteries	220g w/o batteries	>300g incl. rechargeable battery
Size	12x13x5.7cm	12x13x5.7cm	14.8x6.6x4.6cm

\*Clinical studies, PersonalFit Flex compared to PersonalFit data on file (NCT02496429, NCT02492139, NCT02719548, NCT03091985)

You can find more help and advice about breast milk feeding and Medela's pumps and other products at [www.medela.com/flex](http://www.medela.com/flex)